

UTAS Patient Monitoring Solutions

Patient Safety is a Major Priority



Solutions that support patients
and assist clinicians



Why patient safety is so important?

Patient safety is now being recognized as a large and growing global public healthcare challenge.



- The World Health Organization reports that medical errors are the third leading cause of death.¹
- In a recent Philips-sponsored survey emphasize that incomplete patient data, inconsistent care, and alarm fatigue are the main threats to patient.²

Enhanced monitoring solutions can improve patient safety by providing complete patient data and continuous gap-free monitoring.



- UTAS patient monitoring solution is a single monitoring platform for all clinical areas. Move from episodic monitoring to seamless surveillance through our responsive platform concept.
- Modern software and hardware architectures provide flexible integration into hospital workflow.
- UTAS Patient monitors provide real-time vital signs, clinical patient information, and decision support tools at the point-of-care.

To prevent diagnostic errors and guarantee patient safety

- UTAS patient monitors provide advanced functionality and workflow efficiency.
- Innovative multi-parameter alarm analysis reduces false alarms and alarm fatigue.



Protect your investment by choosing state-of-the-art UTAS Patient monitoring solutions

- Unsafe health care can have tragic consequences for individual patients, but its effects reach to the global problem: a lack of focus on patient safety has major financial issues. Usage of low-quality cheap patient monitors resulting medical errors — cost millions!³
- The reliability of the parameters increases the accuracy of the diagnosis, therefore reduces the probability of iatrogenic diseases⁴.
- Enhanced patient monitoring can save lives and money.



High-quality vital signs measuring, decision support tools and personalization of care are becoming a quality standard. UTAS monitoring solutions provide clinical and economic value as well as safety of life.

¹WHO Report by the Director-General : Patient safety Global action on patient safety - 2019

²Patient safety in critical condition : Patient safety survey – 2017 - Phillips/Regina Corso Consulting, 2017

(John H. Eichhorn, MD. Standards for Patient Monitoring During Anesthesia at Harvard Medical School/ John H. Eichhorn, MD; Jeffrey B. Cooper, PhD; ³WHO guidelines for safe surgery 2009 : safe surgery saves lives. Geneva: World Health Organization; 2009, accessed 13 February 2019.

⁴Krishnan NR, Kasthuri AS. Iatrogenic Disorders. Med J Armed Forces India. 2005 Jan;61(1):2-6. doi: 10.1016/S0377-1237(05)80107-8. Epub 2011 Jul 21. PMID: 27407693; PMCID: PMC4923397.

We offer universal monitoring solutions to emphasize individuality

To ensure patient safety, it is important to keep in mind that each patient and his case are unique.



- Therefore, the patient monitoring parameters should be customized according to individual patient's conditions.
- UTAS patient monitors provide extremely flexible modular structure due to **UniPort™ multipurpose connectors**' technology.

Parameters for every clinical situation

- Basic set: ECG, HR, RR, NIBP Smart, Temp, SpO₂ (Masimo SET® or OxiMax™), CO₂ (main- or sidestream*).
- With universal connectors UniPort™ easily attach a wide range of extension modules-in-cable. They will be automatically detected by patient monitor.
- Easy plug and play connection — maximum flexibility and readiness to adjust to new clinical situations.



Extended set of UniPort™ modules:

- | | |
|---------------|--------------|
| • ECG 12-lead | • BIS™ |
| • IBP | • C.O. |
| • Multigas | • ICG |
| • NMT | • EEG / αEEG |

Multipurpose connectors technology enable clinicians to create the most effective configuration of monitoring — **no more and no less as it is necessary for each clinical situation and each patient**



Continuous gap-free monitoring

Seamless monitoring is a complete cycle of **recording**, **transmission** and data **visualization**: from the point of care to patient discharge.



Patient data flow — integrate monitoring data and clinical networked information in one patient record by **UniViewer technology**.


Convenient visualization — view all patient data in the most appropriate way:

- **UniScreen™** extended visualization
- **Bed-to-Bed** solution
- **UNET-S Central station**
- Remote view via **UniViewer**
- Print patients' data on any network printer

Remote view with UniViewer

UniViewer provides remote access to patient data in real-time while you are away from the bedside.

- Access to patient data in the UniViewer network
- Real-time waveforms and parameters
- On smartphones or tablets via browser
- User authorization and data security

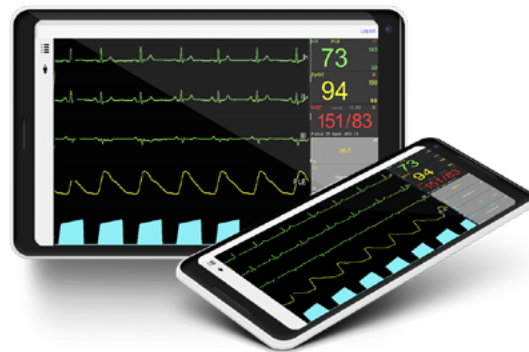

UniScreen™ extended visualization

Customize additional views with UniScreen™ to optimize your workflow.




- Advanced UniScreen™ solution provides the ability to expand and set up visualization effectively on both screens and gives the big picture of patient status.
- UniScreen™ is designed to simplify access to critical patient data you need to make decisions with confidence.



Ambulans,
Transport

Sala
operacyjna

Oddział opieki
po znieczuleniu

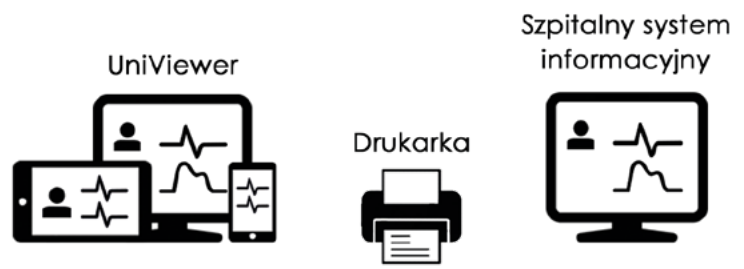

UNET-S Central Station

UNET-S consolidates the monitoring information and provides a quick overview of the patient conditions.


Continues monitoring, interconnection between patient monitors and data recording.

- Integrates up to 32 monitors in wireless or wired connection
- All parameters in real-time displaying
- Continuous parameters recording: trends — up to 720 hours, mini-trends, automatic recording of all parameters when an alarm occurs.
- Patients' database
- Saved data transfer to external medical information systems (HL7)
- Visual and audio alarms



Bed-to-Bed

Get access to patient data from one monitor through any other UM 300-S monitor in the hospital network fast and timely.

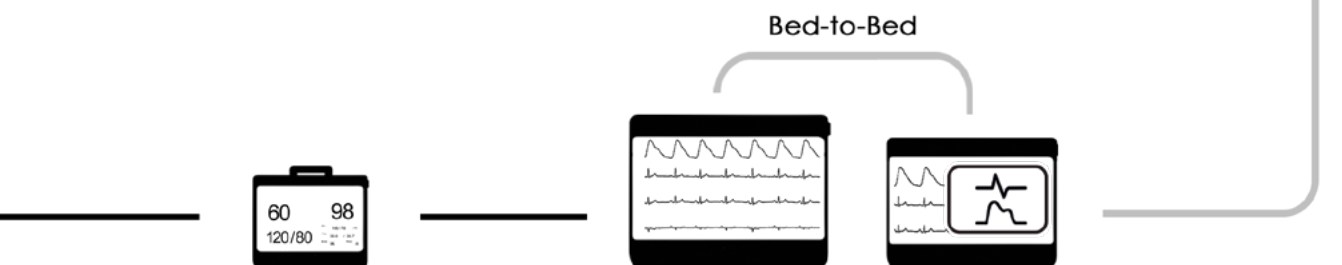

Transport
wewnętrzny



Oddział
intensywnej terapii





Service and support

Be attentive to the state of equipment is important as caring for patients.

We provide a full range of services to maintain superlative state of medical equipment, including remote consultation, diagnostics and software updates for monitors.

To find solutions for even the most difficult situations, we collaborate with leading clinical centers in 30 countries.



We are close to you

Address your question to UTAS and communicate directly with the specialist you need. We maintain consultations in the language of professionals, considering real clinical experience.

At your service are technical experts, product managers, developers, sales & marketing team.

U can save the life®



UTAS Technologies s.r.o.

1, Staviteľská str.,
Bratislava, 83104, Slovakia

tel./fax: +421 220 620 001
tel.: +421 220 620 002
e-mail: info@utastech.eu
www.utastech.eu